

RAPID IMMUNOFLUORESCENT ASSAY FOR TDM* OF BIOLOGICS



Simple

- 1 step to result

Fast

- 10 min for drug measurement
- 12 min for ADAAb measurement

Flexible

- Whole blood (Capillary / Venous)**
- Serum
- Plasma

Reliable

- Excellent correlation with i-Tracker² and USA TRACKER³ reference kits
- Calibrated on NIBSC/WHO international standards**
- Validated on princeps and biosimilars

* Therapeutic Drug Monitoring

** Depending on test: refer to the instruction for use


ez-Track¹

- Time Resolved Fluorescence technology
- Connectable to the LIS
- Friendly user interface
- Compact and lightweight
- Simple and quick test procedure



PROCEDURE

Serum / Plasma

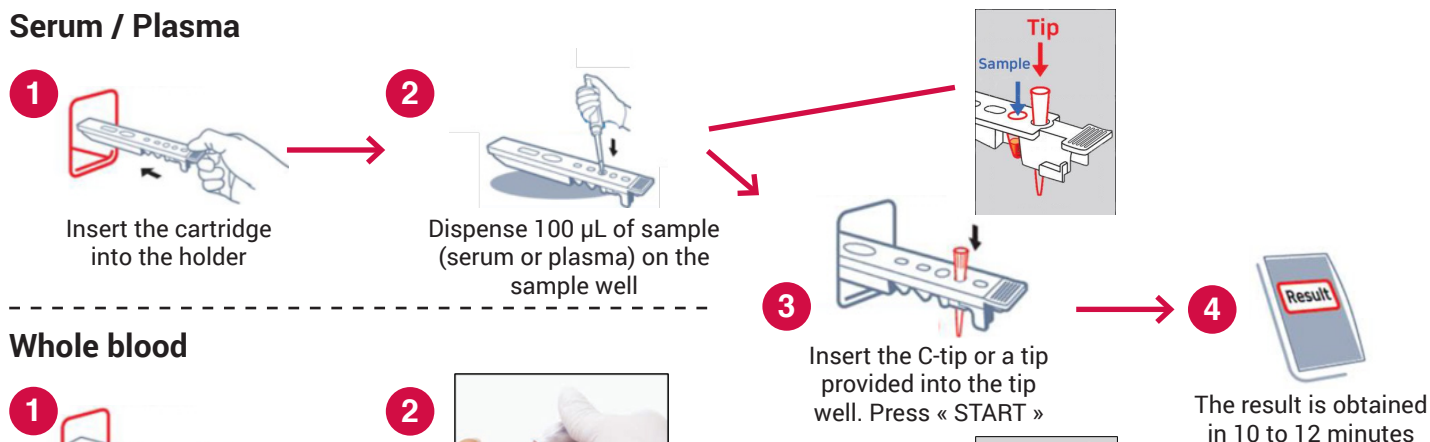

Whole blood



MEASUREMENT RANGE

DRUG		ANTI-DRUG ANTIBODIES	
Infliximab	0.2 - 50 µg/mL	Anti-Infliximab free	4 - 250 AU/mL
		Anti-Infliximab total Ab	8 - 250 AU/mL
Adalimumab	0.2 - 50 µg/mL	Anti-Adalimumab free	3 - 200 AU/mL
Golimumab	0.2 - 30 µg/mL	Anti-Golimumab free	2 - 40 ng/mL
Ustekinumab	0.5 - 20 µg/mL	Anti-Ustekinumab free	2.5 - 100 ng/mL
Vedolizumab	2 - 60 µg/mL	Anti-Vedolizumab free	10 - 250 ng/mL
Etanercept	0.2 - 10 µg/mL	Anti-Etanercept free	10 - 100 ng/mL

PRODUCT REFERENCES

REFERENCE	DESIGNATION	PACKAGING
ET 001	Instrument ez-Track ¹	Instrument and accessories
ETx 002-24	ez-Tracker Drug	24 tests
ETx 003-24	ez-Tracker Anti-Drug Antibodies	24 tests
ETI 003T-24	ez-Tracker Anti-Infliximab Total Ab	24 tests
ETx 002-C / ETx 003-C	ez-Tracker Drug / ADAbs Controls	2 x 1 ml

x = Infliximab / Adalimumab / Golimumab / Vedolizumab / Ustekinumab / Etanercept (Rituximab : in development)

Manufacturer: Boditech (Republic of Korea).
In vitro diagnostic medical devices.
 Carefully read the instructions for use.
 Pictures may differ from actual products.

ez-Track1 - V.06/2024 - EN

BIO SYNEX
 Theradiag

BIOSYNEX Theradiag
 14 rue Amboise Croizat
 77183 Croissy-Beaubourg • France
 Tel. : +33 (0) 1 64 62 10 12
 Fax : +33 (0) 1 64 62 09 66
 E-mail : service.clients@biosynex.com
www.theradiag.com

BIO SYNEX

BIOSYNEX S.A.
 22 Boulevard Sébastien Brant
 67400 Illkirch-Graffenstaden • France
 Tel. : +33 (0) 3 88 78 80 88
 Fax : +33 (0) 3 88 78 76 78
 E-mail : service.clients@biosynex.com
www.biosynex.com